


Channel disable for Bragg adjust.

Several methods.

1: From **Signal Path** menu

If no External modulation signal is connected to J7, select **Int/Ext** button to turn off channel On/Off
RED = OFF

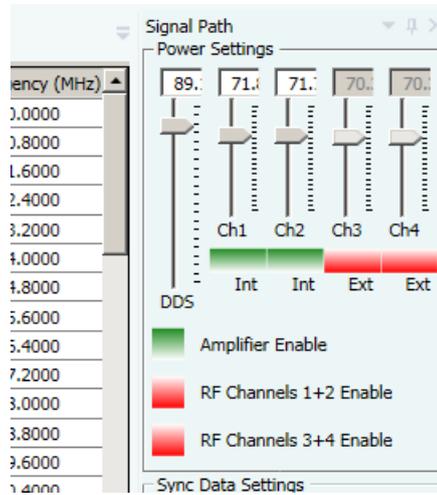


2: From **Calibration** menu

Step 1: IMPORTANT

First reduce the amplitude to ZERO



Step 2: Then select the HOLD button(s) under the channel(s) to turn OFF

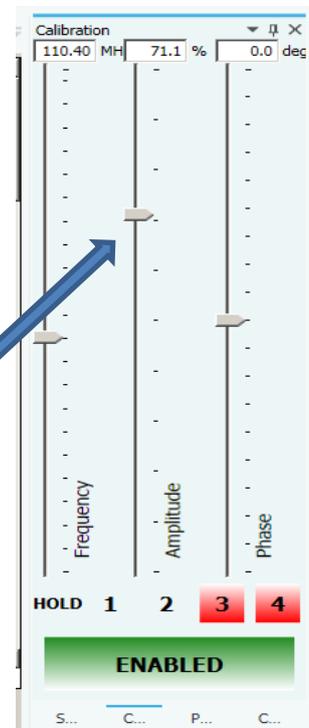

Step 3: Adjust the amplitude for the channels enabled (In this case Ch1 and Ch2)

3: Physical disconnect

Disconnect the coax cables from the IMS4 output(s); J1, J2, J3 and/or J4